



**We recognise
the rhythm!**

© custo·med
EXCELLENCE IN DIAGNOSTICS

custo guard holter with custo wing

Modern Holter with novel patch electrodes

custo guard holter is our modern Holter ECG recorder that combines the benefits of both long-standing experience and continuous development. It excels in high patient comfort, is compact and lightweight, and comes without any disturbing cables or belts. With three channels, motion and pacemaker detection as well as an almost unlimited recording time of several months, the custo guard holter is the perfect solution for an excellent long-term ECG.

custo guard holter: The perfect all-in-one Holter monitor!


**Easy.
Wireless.
Reliable.**





Your benefits:

- 3 channels with two independent leads
- Completely wireless and perfect for hygienic application
- Optimum recording quality with the custo wing electrode
- Fast and stable device communication and data transfer
- Short battery charge times
- Pacemaker detection
- Position and motion detection
- Check for correct application through PC monitoring
- Bluetooth Low Energy for mobile applications
- LED status display for charge status and HR detection
- Protection class IP65: perfect for disinfection
- Very long recording time (several months)
- Running time up to 7 days with only one battery
- Very small and lightweight with only 30 g
- Docking station custo base with one or six spots for device communication and as charging station

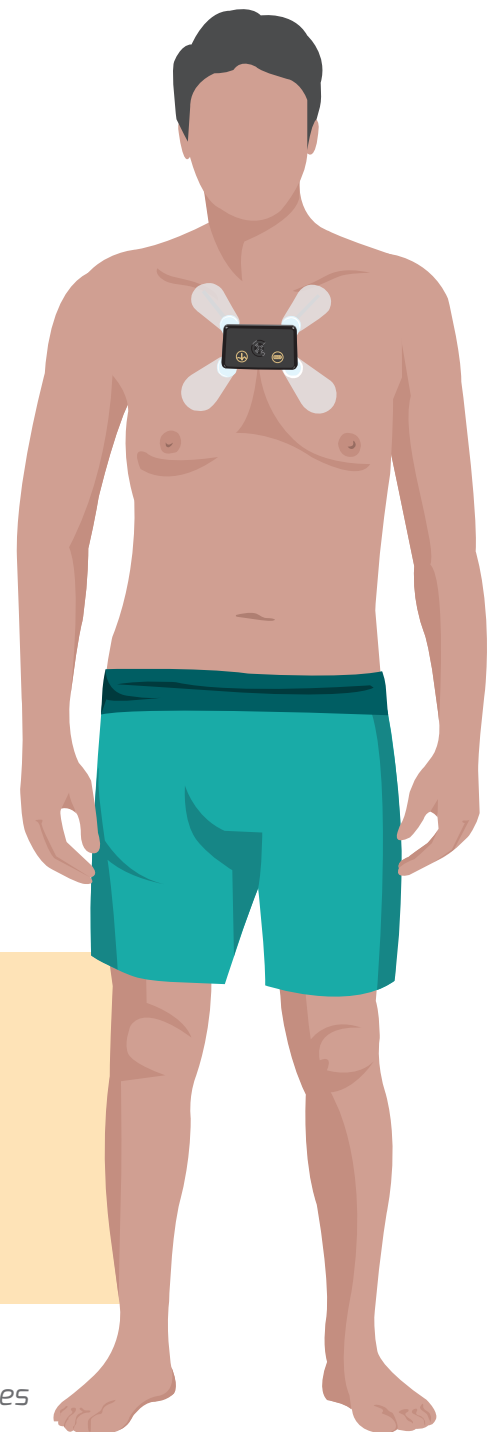

Various electrode applications are possible:

- Special disposable electrode custo wing: Artefact reduction and signal optimisation due to decentral connections and fully covered adhesive surface – material made of breathable nonwoven fabric
- ECG cable guard 3 and ECG cable guard 4
- custo belt 3

Intelligent interface (guard holter viewer) in custo diagnostic

- Display for recording stand-by mode and particular device states
- Various configurations and recording modes are possible

Best in cardiopulmonary diagnostics for more than 40 years!


Since 1982, our brand has been representing innovative system solutions in cardiopulmonary diagnostics. In cooperation with committed and competent sales and service partners, we offer professional solutions taking account of your individual organisational circumstances and IT structures. No matter if your installation is based on a stationary or mobile system, a single-user workstation or multi-site network solution, in a doctor's office or in a hospital – your modular custo med system will be a profitable investment for the future.

Our modular platform custo diagnostic is the intelligent link between clinical engineering and administrative software.



Visit us
on YouTube!

www.customed.de

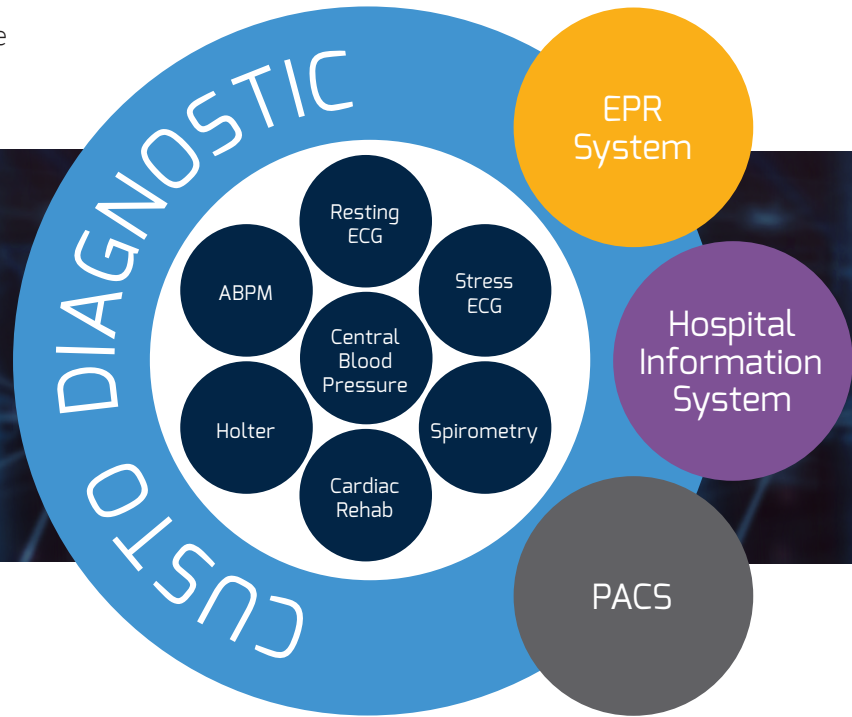

custo med – your local dealer

Contact:
custo med GmbH
Maria-Merian-Str. 6
85521 Ottobrunn
Germany

Phone: +49 89 71098 00
Fax: +49 89 71098 10
Web: www.customed.de
E-mail: info@customed.de

