

NEW

custo diagnostic 5.8

Sport ECG...What is still considered as normal?




© custo·med
EXCELLENCE IN DIAGNOSTICS



NEW: custo diagnostic 5.8 *Sport ECG as an option in the resting ECG*

Functional and structural adaptations of the heart as well as physiological adaptations of the autonomic nervous system can lead to training-related changes in the resting ECG of athletes. These changes in the ECG may cause diagnostic difficulties in such a case.

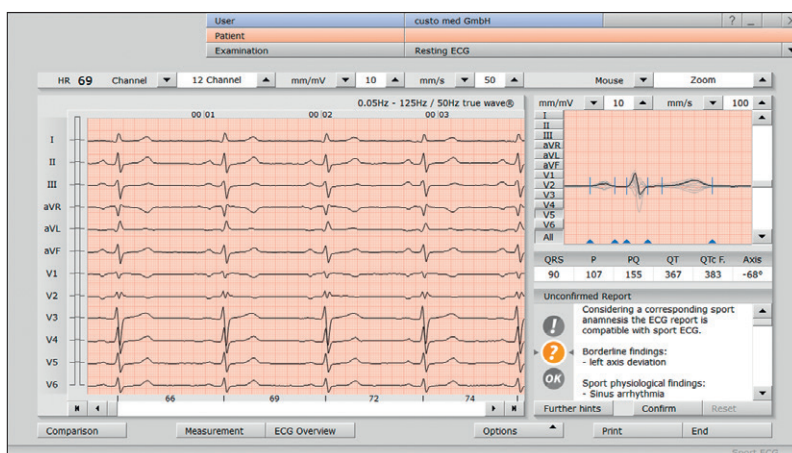


The sport ECG as an option in the resting ECG of the new custo diagnostic 5.8

The question arises: What is still considered as normal? Therefore, the aim of the sport ECG option in custo diagnostic is to differentiate between sport-related physiological and pathological ECG changes by means of a traffic light system. In accordance with the current recommendations of the ESC, normal (training-related) ECG changes are distinguished from borderline and abnormal (not training-related), and thus possibly pathological ECG changes that require further clarification.



With the sport ECG option, the results of the automatic analysis in the new custo diagnostic 5.8 are evaluated differently for competitive athletes than for people who do not fall into this category. The interpretation is based on the “International recommendations for electrocardiographic interpretation in athletes” of the ESC.

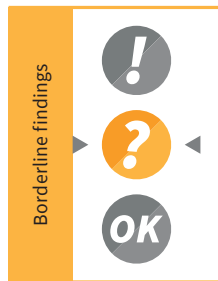

The clear traffic light display is the centrepiece of the sport ECG module.

Resting and stress ECG system custo cardio 300

custo cardio 300 is our 12-channel resting and stress ECG for mobile use. Maximum flexibility is guaranteed due to wireless operation and the possibility of working with one unit at several workstations. Not least, custo cardio 300 is the perfect solution thanks to its modular operability (various patient cables, connection to a suction unit, etc.).



Resting and stress ECG system custo cardio 400/400 accu

custo cardio 400/400 accu sets new standards in functionality, usability and hygiene. Due to the easy-to-change, reusable suction electrodes a hygienic and sustainable application is guaranteed. The version with rechargeable battery is completely self-sufficient and can be wirelessly integrated into your WiFi network, making it flexible to use.



Innovative system solutions for your diagnostics


Since 1982, our brand has been representing innovative system solutions in cardiopulmonary diagnostics. In cooperation with committed and competent sales and service partners, we offer professional solutions taking account of your individual organisational circumstances and IT structures. No matter if your installation is based on a stationary or mobile system, a single-user workstation or multi-site network solution, in a doctor's office or in a hospital – your modular custo med system will be a profitable investment for the future.

Our modular platform custo diagnostic is the intelligent link between clinical engineering and administrative software.



Visit us
on YouTube!

www.customed.de

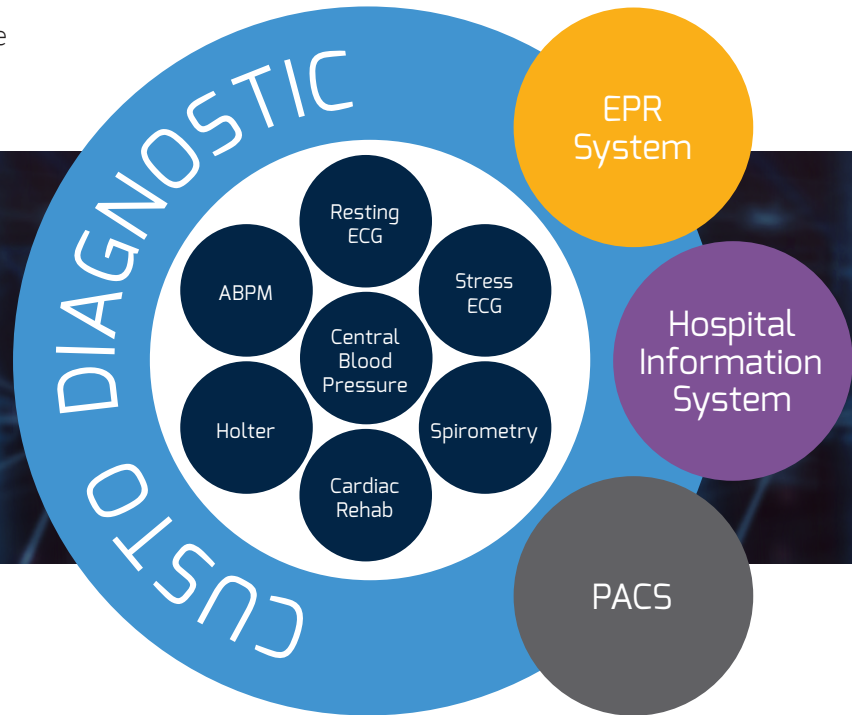

custo med – your local dealer

Contact:
custo med GmbH
Maria-Merian-Str. 6
85521 Ottobrunn
Germany

Phone: +49 89 71098 00
Fax: +49 89 71098 10
Web: www.customed.de
E-mail: info@customed.de